



Agilent CytoGenomics 1.0

Installation and Administration Guide

Research Use Only. Not for Diagnostic Procedures

Installing the Software	2
Flow chart for installation	2
System Requirements	4
Administering the System	5
To manage Users and Roles	5
To manage the Database	6
Getting Help	8
To contact Agilent Technical Support	8
To get help from the Help guides	8
To learn about Agilent products and services	8

The Agilent CytoGenomics software takes little time to install if you prepare properly. And only a few tasks are required of Administrators, although the role of Administrator allows you to perform all the tasks within the program. This guide shows you how to install and administer the program so users can work with the program most effectively.


See the *Product Overview Guide* for an overview of the program and how the tasks within it are organized by roles.




Agilent Technologies

Installing the Software

Flow chart for installation

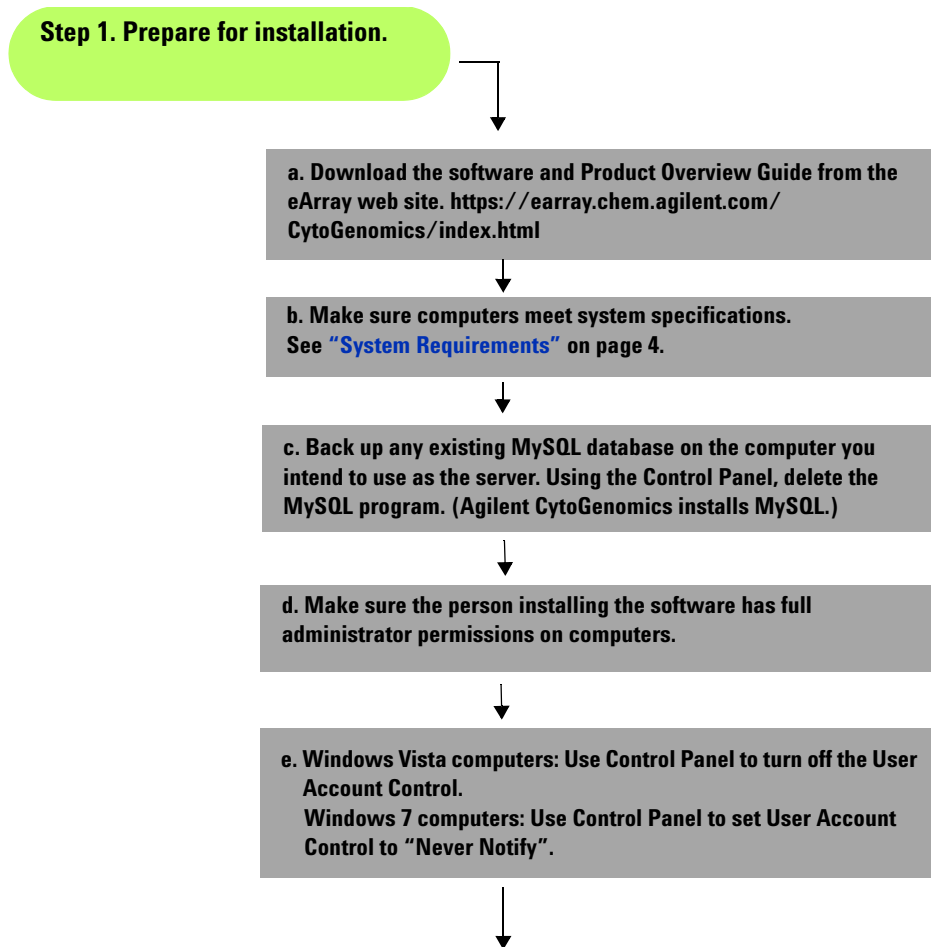

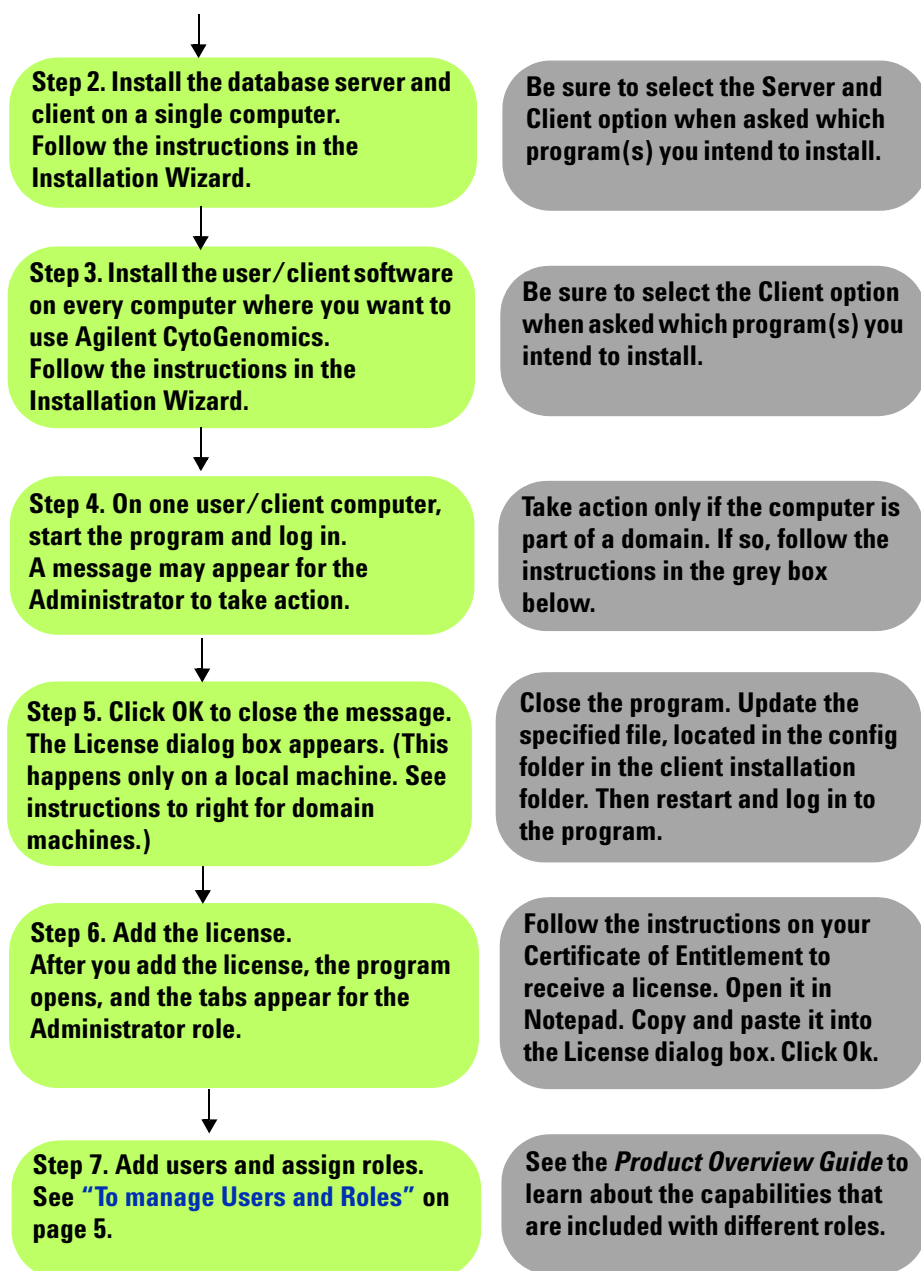

System Requirements

Check to make sure that your computer systems meet the following requirements. All client installations point to the computer on which the server software is installed. See [Table 1](#) for minimum requirements. See [Table 2](#) for recommended requirements.

Table 1 Minimum Requirements for Agilent CytoGenomics computers

Operating system	32-bit Windows (XP-SP2, Vista Business, Windows 7 Professional) If the operating system is not in this list, it is not supported.
Processor	> 2GHz
Available memory (RAM)	4GB
Hard disk space	500GB
Database Software	If MySQL is already installed, you must delete it. Agilent CytoGenomics will not install properly if MySQL is present. Agilent CytoGenomics installs and configures MySQL automatically.
Display Resolution	1280 x768 minimum

Table 2 Recommended Requirements for Agilent CytoGenomics computers (only if not minimum)

Operating system	64-bit Windows (XP-SP2, Vista Business, Windows 7 Professional) systems are supported. If the operating system is not in this list, it is not supported.
Processor	> 3GHz
Available memory (RAM)	8 GB
Hard disk space	Database machine: 500 GB of dual RAID hard drive Client machines: 500 GB

Administering the System

As an Administrator for Agilent CytoGenomics, you have access to everything! That is, you can perform any task available in the program. Your primary responsibility is to add users and their roles to the system, and you also manage the database.

See the *Product Overview Guide* for an explanation of the tasks that each role can perform with the program. Users who are assigned the role of Technician can only run workflows and view the results. They cannot set up workflows, nor can they sign off on results. Those are tasks for users assigned the role of Scientist.

This section shows you how to add users, assign roles and manage the database.

To manage Users and Roles

Only the Administrator role can perform these tasks.

Table 3 Tasks for Managing Users and Roles

To do this	Follow these instructions	Comments
<div>Add new users</div>	<div><div>1 Click Admin, and click the Users button.</div><div>2 At the bottom of the screen, click Add New User. A list of all the users in the domain or on the local machine appears.</div><div>3 Type in the User Name, and click Find.</div><div>4 At the bottom of the Add User Dialog dialog box, click Add User. The user now appears in the list of users with the role of Technician, which has been enabled.</div></div>	<div>A user cannot be deleted but can be disabled.</div> <div><div>a Click Edit for the user you intend to disable.</div><div>b Clear the Enabled checkbox, and click Save.</div></div>
<div>Change roles</div>	<div><div>1 On the line of the user whose role you intend to change, click Edit.</div><div>2 Highlight the role on the right of the panel, then click the < button to remove the role on the right.</div><div>3 Highlight a role on the left of the panel, then click the > button to add the role to the right.</div></div>	

To manage the Database

Even though the Administrator role is the only role permitted to change the location of the database or the common storage folder, any user can restart the MySQL service.

Table 4 Tasks for Managing the Database

To do this	Follow these instructions	Comments
<div>Change database locations and info</div>	<div><div>1</div>In the Admin tab click the Database button.</div> <div><div>2</div>Click Change.</div> <div><div>3</div>Type one or more of these:<ul style="list-style-type: none">Database computer namePort numberCommonStorage pathwayWhen you leave the dialog box, the program accepts the new entries.</div>	<ul style="list-style-type: none">If you change the location of your common storage, you must copy the contents of the previous location to the new location manually. Failure to do so may cause unexpected behavior in the program.The location of common storage must be entered using UNC standard notation. (\\machine2\CommonStorage)
<div>Restart MySQL</div>	<ul style="list-style-type: none">If the MySQL service stops, to restart it, click Restart Cyto MySQL Service on the desktop.	

Getting Help

To contact Agilent Technical Support

Technical support is available by phone and/or e-mail. A variety of useful information is also available on the Agilent Technical Support Web site.

Resource	To find technical support contact information
Agilent Technical Support Web site	<ol style="list-style-type: none">1 Go to http://chem.agilent.com.2 Select a country or area.3 Under Quick Links, select Technical Support.4 Select from the available links to display support information.
Contact Agilent Technical Support by telephone or e-mail (United States and Canada)	Telephone: (800-227-9770) E-mail: informatics_support@agilent.com
Contact Agilent Technical Support by telephone or e-mail (for your country)	<ol style="list-style-type: none">1 Go to http://chem.agilent.com.2 Select Contact Us.3 Under Worldwide Sales and Support Phone Assistance, click to select a country, and then click Go. Complete e-mail and telephone contact information for your country is displayed.

To get help from the Help guides

To find help for all the other tasks Administrators can perform with Agilent CytoGenomics (those that Technicians and Scientists can perform), click **Help** in the upper right corner of the program window.

To learn about Agilent products and services

To display information about the Life Sciences and Chemical Analysis products and services that are available from Agilent, go to www.chem.agilent.com.

www.agilent.com

In this book

This book gives you instructions to install and administer Agilent CytoGenomics.

© Agilent Technologies, Inc. 2010

Revision A0, November 2010



G2602-90002



Agilent Technologies